


n Page No. 9

EcoRI / Bam HI RDS on 1 Hep 5+6  
ex rxn mixed

10µl 1 DNA  
5µl 10x "B"  
1µl EcoRI  
1µl Bam HI  
33µl H<sub>2</sub>O  
50

Inc 37°C ~ 2 hrs → added 10µl 5x dye to each  
2µl each on 0.7% agarose (1x TBE)

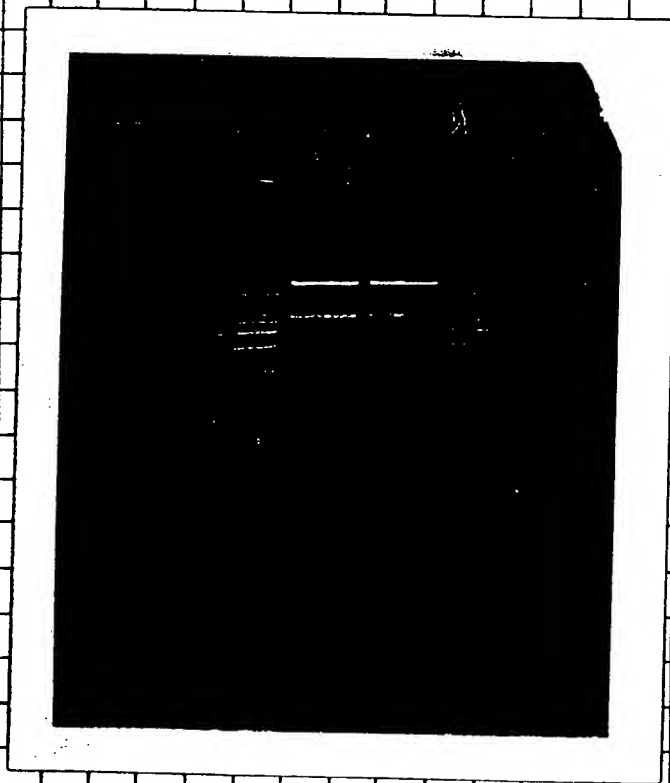

Set up a 1% CMP gel  
for preparative  
electrophoresis on 2/20.  
Apparently 15+6 at  
least have the  
full length coding  
regions for Tyro10.

To Page No. 11

Witnessed &amp; Understood by me,

Date


Invented by

Recorded by

Date

From Page No. 10

Ran remaining 50µl 1546 RD's on 1% agarose  
cut out indicated bands

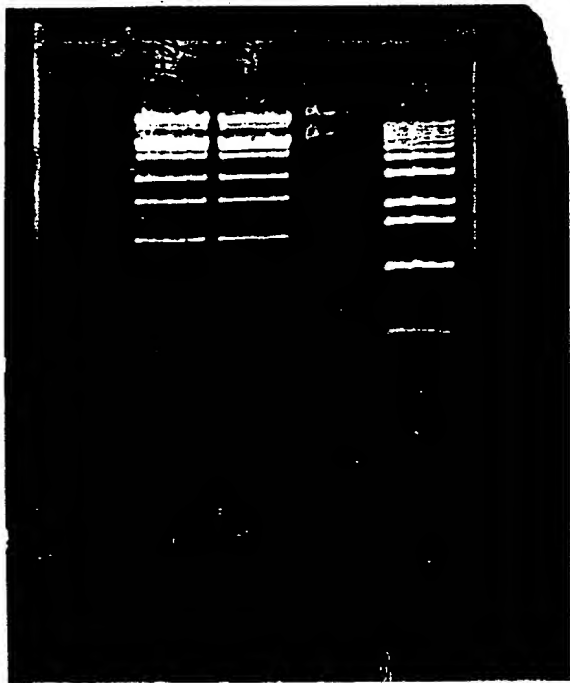

Ran Magic PCR preps  
collected each in 50µl  
combined

Ligation to pRK5  
mixed for each

1µl pRK5 ÷ EcoRI/BamHI / ~~0.001~~  
2µl 10x ligase buffer  
1µl Ligase  
7µl insert (or 1µl)  
8µl H<sub>2</sub>O (or 6µl)

Plus a vector alone ligation  
= 3 total.

Inc 12.5°C O/N.

Witnessed &amp; Understood by me,

Date

Invented by

Recorded by

Date

To Page

Page No. 11

transformed O/N pK5 Ligs to competent  
XL-1 blue cells.

Plated each onto 5 x 100mm LB cant<sup>50</sup> plates  
(cont. plated onto 1 plate)  
Inc all (11) 37°C O/N.

To Page No. 13

ed &amp; Understood by me,

Date

Invented by

Date

Recorded by

WIM Bacon

TITLE

Book No. 175

Exhibit B, pg. 4 of 7

From Page No. 12

RD of SKG / CH<sub>2</sub>CH<sub>3</sub> IgG1 fusion vector

mixed

2 µl DNA (~4 µg)  
 5 µl 10x "B"  
 4 µl H<sub>2</sub>O  
 1 µl EcoRI  
 1 µl BstEII  
 50

Inc 37°C ~ 2 hrs

Ran on 0.7% agarose

Isolated band on N445 paper  
 Eluted in 400 µl 1M NaCl (in TE) 70°C 2.  
 Removed paper.  
 Stored - 70°C o/n.

Ordered PCR primers to PCR ECD of HPI  
 & incorporate a BstEII site @ 3'e

Genentech, Inc.  
Genentech, Inc.

## SYNTHETIC DNA REQUEST

A-6749 ✓

WILL BAEON

2650

10231

1713

LIST SEQUENCE(S) &amp; A NOTE AS SUCH

SIZE &amp; FRAGMENT NAME:

PLEASE INDICATE BY "\*" FRAGMENT(S) TO BE CLONED

1x 24mer  
 1x 35mer

① Tyro P13 (24mer)

⑤' TGG·GAG·GAG·GAG·CCC·ATG·CGC·CAC·3' "

② Tyro P14 (35mer)

⑤' GTA·CAG·TNA·CCG·GCG·GTC·GAG·CTC·  
·CCC·TCG·GAC·TT·3' "

FRAGMENT USE: ☐ PRIMER ☐ PROBE ☒ PCR ☐ GENE CONSTRUCTION ☐ UNKERVADAPTOR  
☐ MUTAGENESIS ☐ OTHER (SPECIFY):

SPECIAL REQUESTS

WILL BAEON

DATE

APPROVED

Parker

DATE

Witnessed &amp; Understood by me,

Date

Invented by

Date

Recorded by

WILL BAEON

To Page N

om Page No. 13

Checked PRK5/HPTK6 (tyro10) transformation plates  
Got 39 transformants

Started 39 x 5ml 1B cult<sup>50</sup> (+ master plate)

Inc all 37°C O/N.

To Page No. 1

Witnessed & Understood by me,

Date

Invented by

11/7/13 Balman

Date

11/7/13

TITLE \_\_\_\_\_

Project No. \_\_\_\_\_  
Book No. 175

Exhibit B, pg. 6 of 7

From Page No. 14

Did Magic MP's on 12 of Q/N. PRK5/HPK6 MP's  
(Stored remaining 27 cultures + master @ 9)

Collected each in 50  $\mu$ l TE

RD's  
Per rxn      2  $\mu$ l MP DNA  
              2  $\mu$ l 10x "B"  
              15  $\mu$ l H<sub>2</sub>O  
              0.5  $\mu$ l EcoRI  
              0.5  $\mu$ l BamHI  
              20

Inc 37°C ~ 1.5 hrs  $\rightarrow$  added 4  $\mu$ l dye to each

Ran ~ 15  $\mu$ l each on 0.7% agarose (1x TBE)



All positives.

Started 1x 500  $\mu$ l ZYT carb<sup>50</sup> Midi  
prep on #1

Inc 37°C O/N.

To Page \_\_\_\_\_

Witnessed &amp; Understood by me, \_\_\_\_\_

Date \_\_\_\_\_

Invented by \_\_\_\_\_

H. H. Benson

Date \_\_\_\_\_

m Page No. 15

torced 1.2ml O/N HPTK6/pRK5 + 400µl 50% glycerol  
E70°C

Did Magic Maxi prep on remainder of culture

EtoH ppt'd DNA after elution from resin column

Stored ppt'd pellet in -20°C 70% EtoH. O/N

To Page No. 17

essed &amp; Understood by me,

Date

Invented by

Recorded by

Date

WIM Bacon